

GALLAY ACCESSORIES

FOR LABORATORY GLASSWARE WASHERS

Basket And Rack Configurations.

Speak to our sales consultants about other configurations to suit your needs.



Two Spraying Levels

Standard lower basket and upper basket with sprayer. In this configuration there are two washing levels with rotating sprayers which can hold all the supports for the various glassware items (including beakers, test tubes, Erlenmeyer flasks and other glassware).



One Injection Level

The hollow jet rack (cat no. LM40) can be considered the universal hollow jet rack because it is ideal for washing narrow-necked glassware of different sizes. It has a capacity of 40 items and can wash both large and small flasks or bottles.



Two Levels- Spraying/Injection

Standard lower basket & top hollow jet rack. In this configuration there are 2 washing levels – a lower level with rotary spraying designed to hold various racks & baskets such as clip or spring peg racks & an upper level 40 position Jet Rack. This configuration allows narrow-necked glassware to be washed at the same time as other glassware.



Two Injections Levels

A rack consisting of two injection levels with a capacity of 68 items. This version delivers optimum injection washing capacity.

 **smeg** INSTRUMENTS for **GALLAY**

FOR ALL YOUR MACHINES AND ACCESSORIES CHOOSE US.

 **GALLAY**
MEDICAL & SCIENTIFIC PTY. LTD.

GALLERY BASKETS AND RACKS...FOR ALL APPLICATIONS

JET RACKS



Hollow Jet Rack LM40

Universal jet rack for washing narrow-necked glassware, round bottom flasks and measuring cylinders. Eight different sizes of interchangeable nozzles. Maximum capacity 40 items.

LM40-DS: Rack with Fan forced Drying.



2-tier Hollow Jet Rack LM80

Two-level rack for washing narrow-necked glassware by internal injection. Capacity 68 items, for glassware with a maximum height of 225mm. Standard equipment includes 140mm height nozzles.

LM80-DS: Rack with Fan forced Drying.



Upper Hollow Jet Rack LM40S

Upper level rack for washing narrow-necked glassware by internal injection. It is positioned on the top level of the machine. Capacity 40 items, for glassware with a maximum height of 225mm. Standard equipment includes 140mm height nozzles.

LM40S-DS: Rack with Fan forced Drying.



Injection Jet Rack for Bottles LB161

Injection Jet Rack for bottles for Micro or Pharma use, available in six versions.

LB161-DS: Rack with Fan forced Drying.



Flask/Pipette Rack LPM20/20

Injection washing rack for both pipettes and flasks. Features 20 positions for pipettes with height measurements of up to 550mm and 20 positions for narrow-necked glassware with a maximum height of 490mm.

LPM20/20-DS: Rack with Fan forced Drying.



Pipette rack LPV40

Rack for washing volumetric pipettes measuring up to 550mm in length. Features 20 positions for 550mm long pipettes and 20 positions for 450-490mm long pipettes.

LPV40-DS: Rack with Fan forced Drying.



Pipette Washer Disk LPT100

Accessory to wash 100 pipettes with 1 to 20 mL capacities and height up to 450mm. The pipette disk must be fitted onto a rack for quick loading and unloading access.

LPT100-DS: Rack with Fan forced Drying.



2-tier Injection Rack for Test Tubes KP200

Two level trolley for washing tapering and straight tubes by internal injection. Maximum 200 tubes. Single level version also available.

KP200-DS: Rack with Fan forced Drying.

1-tier Injection Rack for 100 tubes KP100



Flask Washing Rack LM20

Universal flask washing rack with 20 nozzles on half the trolley. The other half is free for mixed load and can contain 1/2 basket accessories.


LM20-DS: Rack with Fan forced Drying.

BASKETS / RACKS



Upper Basket with sprayer CS1

Fits upper level of washer to support racks of various configurations. Features Spray arm for upper level washing.

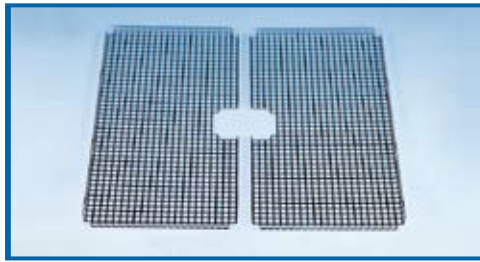


Half Base Grids 7250184

Base Grids - upper with hole 7250182

For Upper Baskets. Fits into upper basket to assist with support of medium to large sized glassware.

Base Grids-Lower 7250183

For Lower Baskets. Fits into lower basket to assist with support of medium to large sized glassware.

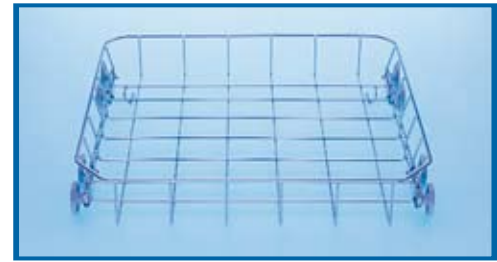

Lower Basket CS2

Rolls into lower position in washer to support baskets or racks of various configurations.



Petri Dish Holders

- PD70S-** Capacity 40 pieces, dia. 50-70mm positioned on upper basket.
- PD100S-** 38 pieces, 70-120mm positioned on upper basket.
- PD70I-** 42 pieces, 50-70mm positioned on lower basket.
- PD100I-** 38 pieces, 70-120mm positioned on upper basket.



Quarter Basket 150mmH 7251550a

Quarter Basket 200mmH 7251552

Dividers 7251550b

Adjustable Lids 7251550c

Lids 7251550d

Please note: This picture displays a set of four (not sold as whole unit)



Spring Peg Rack SB28

Universal Erlenmeyer flask/beaker accessory recommended for capacities up to 1000mL.



22 Clip Rack 7250348

30 Clip rack (not illustrated) 7250349

Clip Racks are used for conical flasks, beakers, measuring cylinders, etc.



Basket for Test Tubes 7250906

Basket (185x108x205mm) with height adjustable lid for test tubes.



Square Basket

205x205x160H mm 7250003

205x205x70H mm 7250004

Lid for above baskets 7250008



Rectangular Basket

410x190x120Hmm 7250006

410x190x70Hmm 7250005

Lids to suit above baskets 7250009



Square Baskets 100x100x100mm

For Test Tubes (no lid) 7251229

For Screw Caps with lid 7250902



Rectangular Basket 7251228a

400x100x100mm

Lids to suit above baskets 7251228b

Adjustable lid to suit above basket 7251228c

GLASSWARE SUPPORT / ACCESSORY SELECTION TABLE

For use in Smeg Laboratory Glassware Washers

Glassware	Spring Peg Racks	Clip Racks	Hollow Jet Racks	Injection Bottle Racks	Pipette Racks	Base Grids	Quarter Baskets	Square Baskets	Rectangular Baskets	Test Tube Injection	Petri Dish Rack	Custom Made Racks
Erlenmeyer Flasks	●	●	●									
Beakers (large)						●		●	●			
Beakers (small)	●	●						●	●			
Measuring Cylinders	●		●							●		●
Volumetric Flasks			●									
Bottles (large)				●								●
Bottles (small)	●		●					●	●			●
Separating Funnels			●									
Funnels		●				●		●	●			
Round Bottom Flasks			●									●
Narrow Necked Items			●	●								●
Tubes and Vials							●	●	●	●		●
Stainless Steel Bowls											●	●
Stoppers and Lids							●	●	●			
Stainless Steel Instruments								●	●			●
Butyrometers			●									●
Pipettes					●							●
Petri Dishes											●	
Watchglass Covers											●	
Quickfit Glassware	●		●					●	●			●

GALLAY CLEANING AGENTS:

AUTOMATIC WASHING MACHINE CLEANING DISINFECTING AGENTS AND RINSE AIDS FOR LABORATORIES AND HOSPITALS



Gallyay Medical & Scientific P/L cleaning agents are especially developed for use in automatic washing machines and efficiently solve problems associated with washing laboratory glassware and hospital instruments. They are readily soluble and emulsify easily to ensure complete elimination of the strongest residues and lime soda. These products may be used with complete confidence on the most delicate laboratory glassware and instruments.

Gallyay Medical & Scientific P/L cleaning products satisfy modern standards of waste control and are all biodegradable.



HEAD OFFICE MELBOURNE
6 Florence Street, Burwood, VIC 3125
Phone: (03) 9888 8422
Facsimile: (03) 9888 8355

BRISBANE OFFICE
Unit 3, 10-12 Hook Street, Capalaba, QLD, 4152
Phone: (07) 3245 3388
Facsimile: (07) 3245 4388

SYDNEY OFFICE
Unit 3/4 Duck Street, Auburn, NSW 2144
Phone: (02) 9647 2002
Facsimile: (02) 9647 2052

Email: sales@gallyay.com.au Web: www.gallyay.com.au